

## Savina Technical Service Bulletin # 3

**Re:** Improvement of EMC

**Date:** December 13, 2001

**Reference Doc:** Immediate report SM 100000000816.


**Reason:** Dräger Medical has released a modification to the Savina Housing P/N 8414280 to improve EMC to 20V/m.

**Devices affected:** Replace the front panel frame on the following 15 devices:  
ARRJ-0013 through -0022  
ARRK-0091, ARRK -0092, ARRK -0095, ARRK -0096, ARRK -0097

**When:** This action is to be completed by end of March 31, 2002

**Reporting:** Report modified Savinas with S/N, labor and traveling time to Mike Kelkart (kelhartm@draegermed.com).

**Additional Info 1:** The new version of the front panel frame has 3 EMC gaskets glued on unpainted areas.

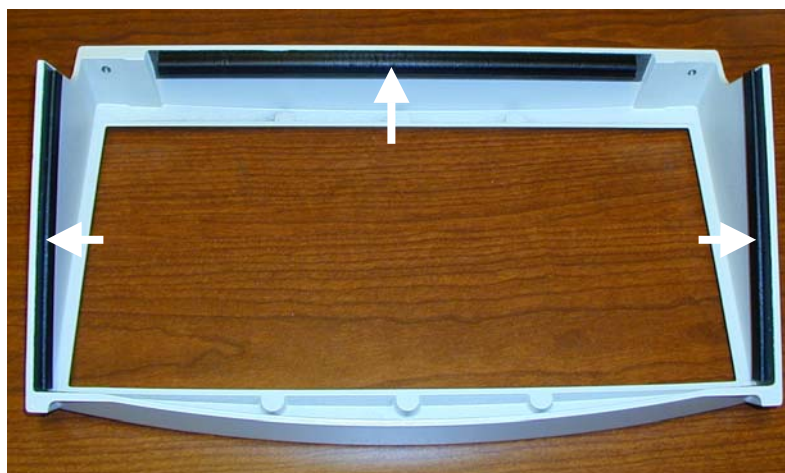

**Additional Info 2:** Replacement time less than 5 minutes.

**Cost:** Free of charge

**Ordering Info:** Front panel frame P/N 8414279  
Note: For this action this part comes with 3 EMC gaskets P/N 8414722 assembled.

**Distribution:** DrägerService Personnel and authorized Service Organizations for CCS products.

If you have any questions, please contact Technical Support by phone at 1-800-543-5047 or by fax at 1-215-721-5789.

Dräger Medical, Inc.  
Technical Product Manager